



Orthopaedics

Allen® Lift-Assist® Beach Chair #A-91500



Lightweight With Lift-Assist™ Technology To Make Positioning Easy

The Lift-Assist Beach Chair is designed for use with orthopaedic procedures using the Fowler's or beach chair position. The chair offers unobstructed access and the Lift-Assist mechanism makes repositioning the patient quick and easy. The Lateral Slide feature enables the surgeon to move the patient closer to the edge of the OR table for increased surgical site access. The Lift-Assist Beach Chair weighs under 30 lbs (13.6 kg) without the back support panels installed.

- Attaches to the rails of standard OR tables
- 500 lbs (227 kg) patient weight capacity
- Lateral adjustment - improved surgical site access
- Lightweight - for single person set-up and removal
- Rail width 20" - 26" (50.8 cm - 66 cm) - fits a wider range of OR tables
- Spring Lock Quick Connect - easy installation and removal of head positioners
- Rail clamp indicators - display when the Beach Chair is properly installed on the OR table rails
- Small rail clamp footprint - frees up accessory rail space needed for different procedures
- No gas spring cables to repair or break
- Removable back panels increases access to surgical site
- Low-cut backboard - excellent maneuverability during total shoulder procedures
- Optional head positioner lock helps support patient's neck



Lift-Assist Pistons Make Positioning Patients Easy



Adjust The Patient Using Lift-Assist™ Handle



Chair Lies Flat For Clear Anesthesia Access

Unobstructed Posterior & Lateral Access



What's Included



#A-91501
Allen® Lift-Assist®
Beach Chair
1 Included



#A-90040
Reusable Pad Set
1 Included



#A-91502
Intubation Plate
and Pad
1 Included



#A-91503 (US)*
Lift-Assist Beach
Chair Clamp
2 Included

What's Required

Counter Traction



#A-90015
Counter Traction Strap
1 Required

+



#A-40022 (US)**
Easy Lock® Socket
1 Required

Headrest & Disposables



#A-90026
Universal™ Head
Positioner
1 Required

+



#A-90023
Universal Head
Restraint Disposables
Case of 10

Arm Support *Choose One Arm Support*



#A-90002
LPS Arm Support
1 Required

OR



#O-ABMT
Multi-Task Armboard
1 Required

+



#O-AMA97X1.0S
Clamshell Armboard
Pad 1" (2.5 cm)
1 Required

+



#A-40022 (US)**
Easy Lock Socket
1 Required

What's Optional



#A-90100 (US)**
Beach Chair Cart



#A-92000
Allen Arm Positioner



#P2000
Thermedx® FluidSmart™
Fluid Management System



#O-LPA
Allen Easy Irrigation Tower

*Specify Rail Type

US #A-91503
EU #A-91503EU
UK #A-91503UK

Japan #A-91503JP
Denyer #A-91503DEN

**Specify Rail Type: US, EU, UK, JP, DEN. See product page for international part numbers.



Universal™ Head Positioner



Adjusts Intraoperatively



Easy To Use Disposable



Single Knob
Releases All Joints

Universal™ Head Positioner #A-90026

Easy To Use And Protective Head Positioning For The Beach Chair Position On The Operating Room Table

The Allen® Universal Head Positioner supports the patient's head and neck providing adjustability and ease of positioning. Designed for simple set-up and use. It supports the patient's unique anatomy securely and comfortably and offers more versatility than competitive head supports.* In addition, the disposable straps which support the head and neck are comprised of wide soft protective padding. Straps are designed to distribute force to help reduce the risk of pressure injuries.

Ideal for use with most upper extremity procedures in the beach chair position.

- Multiple straps for secure control of patient's head and neck
- Single knob releases all joints
- Easy intraoperative adjustments
- Adjustable to help accommodate the patient's head
- Disposable wide soft straps minimize pressure on skin
- Lightweight and portable
- Simple set-up, easy to use
- Secure head and neck with ability to lock the desired position
- Compatible with all Allen Beach Chairs
- Device Weight (without disposable): 3.4 lbs (1.5 kg)
- Patient Weight Capacity: 500 lbs (227 kg)

What's Included



#A-90026
Universal Head Positioner
1 Included

What's Required



#A-90023
Universal Head
Restraint Disposable
Case of 10

What's Optional



#A-91500
Lift-Assist® Beach Chair



#A-90100 (US)**
Beach Chair Cart



#A-90025
Replacement Straps

*Based on a feature comparison to Linvatec®, McConnell and Mizuho OSI beach chair head rests. (January 22, 2010)

**Specify Rail Type: US, EU, UK, JP, DEN. See product page for international part numbers.



Unobstructed Access



Multiple Positions



External Rotation & Abduction

Simple &
Lightweight



Allen® Arm Positioner #A-92000

Spend Your Time On The Procedure, Not The Positioning

The Allen Arm Positioner delivers easy arm positioning for shoulder surgery in the Beach Chair position, as well as elbow and clavicle procedures. By indexing at the elbow, the device is able to gain complete positioning control of the shoulder joint. The revolutionary patent pending design allows the entire device to be steam sterilized and provides an unobstructed access to the surgical site. The device also has the capability of achieving a semi-locked state allowing the surgeon to simply move the arm where he wants it without having to adjust the tension.


Ideal for use with all open and arthroscopic upper limb procedures.

- No foot pedal, no hoses, no hydraulics, no power cords
- The only fully autoclavable arm positioner on the market
- Single knob releases all joints
- Lightweight device weighs only 6 lbs (2.7 kg). No cart required.
- Easy, fluid positioning for enhanced site access
- Connects over the drape, in the sterile field
- Compact and portable for easy storage
- Index at the elbow for complete control
- Quick and easy set-up
- Unobstructed access to the surgical site (no bars or foot pedal obstructions)
- Patient Weight Capacity: 500 lbs (227 kg)

What's Included



#A-92000
Allen Arm Positioner
1 Included

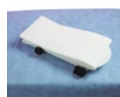

#A-92001
Allen Arm Positioner
Sterile Disposable
Case of 10



#A-92002
Autoclavable Clamp
1 Required

What's Optional



#A-91500
Lift-Assist® Beach Chair



#A-90100 (US)*
Beach Chair Cart




#A-92003
Replacement Arm Support

*Specify Rail Type: US, EU, UK, JP, DEN. See product page for international part numbers.



Multi-Task Armboard #O-ABMT

Provides arm support for lateral, chair, prone and neuro park bench positions. Lockable ball socket combined with 6" (15 cm) of horizontal travel permits a wide range of positioning flexibility. Patient Weight Capacity: 500 lbs (227 kg).

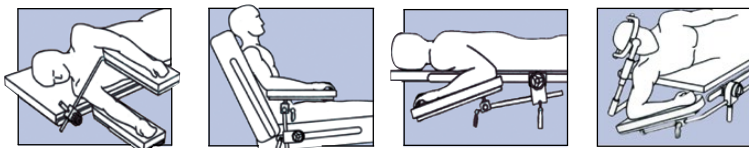

What's Required *Choose One Pad*



#A-40022 (US)*
Easy Lock* Socket
1 Required



Clamshell 1" (2.5 cm)
Armboard Pad
#O-AMA97X1.05
1 Required

OR



Standard 1" (2.5 cm)
Armboard Pad
#O-AMA12X1.05
1 Required

*Specify Rail Type: US, EU, UK, JP, DEN. See product page for international part numbers.



The Wrap-Around Pad Provides
Protection & Stability

LPS Arm Support #A-90002

Provides comfortable and convenient placement of non-operative arm. Easy access to patient IV and enhanced stabilization of shoulder. Patient Weight Capacity: 500 lbs (227 kg).

What's Included



#A-90002
LPS Arm Support
1 Included



#A-AMA74
LPS Arm Support Pad
1 Included

What's Required



#A-40022 (US)*
Easy Lock* Socket
1 Required

*Specify Rail Type: US, EU, UK, JP, DEN. See product page for international part numbers.



FlexiForm™ Head Rest Pad Disposables

#A-90019

This disposable pad creates a more comfortable and protective environment for your patient. Required for use with the FlexiForm Head Rest.
(Sold In Cases Of 10)



Universal Head Restraint Disposables

#A-90023

The Universal Head Restraint Disposable is designed for use with the Allen® Universal Head Positioner and Allen beach chairs during shoulder procedures, specifically in the Fowler's position. The disposable supports the patient's head and jaw helping align the neck in a natural position. (Sold In Cases Of 10)



Allen Arm Positioner Sterile Disposables

#A-92001

The Allen Arm Positioner Sterile Disposable is used exclusively with the Allen Arm Positioner. (Sold In Cases Of 10)



Counter-Traction Brace #A-90001

Applied to operative side, the Counter-Traction Brace stabilizes patient for shoulder manipulation. For use with the Lift-Assist® Beach Chair.

What's Required



#A-40022 (US)*
Easy Lock® Socket
1 Required

*Specify Rail Type: US, EU, UK, JP, DEN. See product page for international part numbers.



Counter-Traction Strap #A-90015

The Counter-Traction Strap supports the patient while they are in Fowler's position. For use with the Lift-Assist Beach Chair.

What's Required



#A-40022 (US)*
Easy Lock® Socket
1 Required

*Specify Rail Type: US, EU, UK, JP, DEN. See product page for international part numbers.



Beach Chair Reusable Pad Set #A-90020

Pad Set For Use With The Allen® Beach Chairs #A-91000 & #A-90000

Combines high-density foam with a durable 4-way stretch cover for patient comfort.



Lift-Assist® Beach Chair Reusable Pad Set CE

#A-90040

Pad Set For Use With The Lift-Assist Beach Chair #A-91501

This pad set combines durable high-density foam with a 4-way stretch cover for patient comfort.



Allen Beach Chair Cart #A-90100

The perfect companion, the Allen Beach Chair Cart conveniently stores the Allen Beach Chair and all of its accessories. The Beach Chair rests securely on sturdy rails at user height to help reduce potential back strain. Optional accessory hooks allow quick and easy storage for larger accessories. Deep bin affords plenty of extra storing space. Weighted base makes cart tip-resistant and reduces potential vibration during transport. Polymer base is dent-resistant. Heavy-duty casters roll smoothly. Two casters lock for safer and easier loading, unloading and storage.

What's Optional



Specify Rail Type
 #A-30011 (US)
 #A-30011E (EU)
 #A-30011UK (UK)
 #A-30011J (Japan)

The Accessory Hook allows for faster product storage on our Accessory Cart. Coated prongs help protect stored product from developing scratches. Easy Lock® design makes attachment and detachment easy.




#A-30012
Single Hook



#A-30013
Double Hook

Specify Rail Type

US	#A-90100	Denyer	#A-90100D
EU	#A-90100E	Japan	#A-90100J
UK	#A-90100UK		

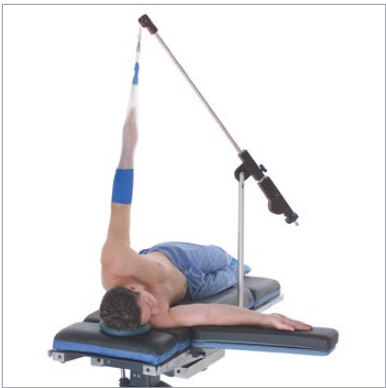

Weightless Shoulder Suspension System

#R-644-36

The Weightless Shoulder Suspension System allows for approximate traction settings up to 40 lbs (18.2 kg) with a simple turn of the tension adjustment knob. Boom arm adjusts from 0° to 90°. The boom head swivels so that its pulley always aligns properly with the cord. The system folds down for easy storage. Includes one sterile Arm Trap Disposable.

What's Required



#A-40018 (US)*
Simple Clamp
1 Required



#A-21200
Sterile Arm Trap Disposables
Case of 6

What's Optional



#A-30010 (US)*
Accessory Cart



#O-UPBL3S
Lateral Brace
3 Piece Set




#A-60003, #A-60002
Allen® Hug-U-Vac®
Lateral Positioner
Standard or Large



#O-LPA
Allen Easy
Irrigation Tower

*Specify Rail Type: US, EU, UK, JP, DEN. See product page for international part numbers.

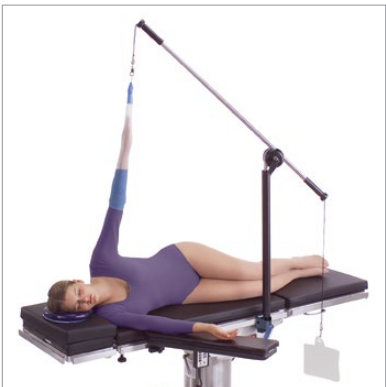

Weights Not Included

Arm Traction Device #A-20500

Accommodates intraoperative traction up to 40 lbs (18.2 kg) for shoulder procedures. A mechanical safety stop helps prevent patient injury in the event the support arm is inadvertently released. Save time and cost by using it as an arm support for both prepping and traction.

What's Required



#A-40018 (US)*
Simple Clamp
1 Required



#A-21200
Sterile Arm Trap Disposables
Case of 6

What's Optional



#A-30010 (US)*
Accessory Cart



#O-UPBL3S
Lateral Brace
3 Piece Set



#A-60003, #A-60002
Allen® Hug-U-Vac
Lateral Positioner
Standard or Large



#O-LPA
Allen Easy
Irrigation Tower

*Specify Rail Type: US, EU, UK, JP, DEN. See product page for international part numbers.



Lateral Braces #O-UPBL3S

Extended posts adjust for use on a variety of patients. Stabilizes patients on standard OR tables. For total hip and other procedures requiring optimum lateral stability. Patient Weight Capacity: 400 lbs (181 kg).

- Braces are fully adjustable
- Soft pads provide soft cushioning of bony prominences
- Swivel heads conform to anatomical contours of patient and can be securely locked for precise positioning

What's Included



#O-UPBLR
Single Rectangular Support
2 Included



#O-UPBLC
Single Circular Support
1 Included

What's Required



#A-40040 (US)*
Easy Lock® Blade Clamp
3 Required

What's Optional



#A-30010 (US)*
Accessory Cart

*Specify Rail Type: US, EU, UK, JP, DEN. See product page for international part numbers.



Sterile Arm Trap Disposables #A-21200



Made of soft plastic mesh, the Allen® Sterile Arm Trap comfortably holds the arm. Arm trap is affixed to patient with CoFlex® NL Self-Adherent Wrap (roll included). Rope loop at the top of the trap allows for easy attachment to many types of traction and suspension devices. (Sold In Cases Of 6)



Integral Straps



Safety Vacuum Valve



Waterproof Disposable Slipcover



Allen® Hug-U-Vac® Lateral Positioners

#A-60003, #A-60002

Safe, Comfortable And Secure Lateral Patient Positioning

The Allen Hug-U-Vac Lateral Positioner safely and comfortably holds patients in lateral and supine positions. The Allen Hug-U-Vac Lateral Positioner's integral quick-release straps enable it to attach securely to the rails of standard OR tables.

The Allen Hug-U-Vac Positioner has a secure valve that eliminates the need for continuous suction.

Ideal for use with Urology Procedures, Orthopaedic Procedures, Cardiothoracic Procedures, Neurological Procedures, Thoracotomy Procedures, Plastic Surgery and General Surgery. Can be used for supine positioning in a variety of procedures.

- Enables secure lateral and supine positioning in a variety of procedures
- Integrated axillary pillow
- Lightweight for easy transport, set-up and storage
- Valve location designed for quick access to anesthesia
- No continuous suction required
- Durable vacuum valve with on/off switch and lock
- Disposable slipcovers are waterproof for easy clean-up

#A-60003 Standard

Patient Weight Capacity: 500 lbs (227 kg)
Device Dimensions: 36" x 32" (91.4 cm x 81.3 cm)

#A-60002 Large

Patient Weight Capacity: 800 lbs (363 kg)
Device Dimensions: 48" x 35" (121.9 cm x 88.9 cm)

What's Included *Choose One*



#A-60003
Standard Allen Hug-U-Vac
Lateral Positioner
1 Included



#A-60002
Large Allen Hug-U-Vac
Lateral Positioner
1 Included



#A-60103
Standard Allen Hug-U-Vac
Lateral Disposable Covers
Case of 12



#A-60102
Large Allen Hug-U-Vac Lateral
Disposable Covers
Case of 12

What's Required *Choose One*

What's Optional



#O-RSX2
Allen Large
Security Strap



#C-008001
Replacement
C-Clip Hug-U-Vac



Active Heel™ Technology
Shown Without Boot Shell



One Boot Accommodates Many
Foot Sizes & Left Or Right Foot



Easy To Use Buckles



Allen® Traction Boot

Active Heel™ Technology Provides Superb Security To Help Prevent Foot Slip

Positioning the leg and foot properly during trauma and hip fracture cases is critical to a less stressful surgery and successful outcome. Many hospitals are still using the original boots that may be decades old and require copious amount of tape to help secure the foot. To mitigate these challenges, Allen has introduced the Allen Traction Boot with Active Heel Technology to help prevent foot slip while providing superb security.

With its modular design, you can upgrade your existing trauma tables to the latest technology that has benefits for the patient, surgeon and OR staff.

Active Heel technology utilizes two cushioning pads on the anterior and posterior portions of the leg to provide a secure grip on the entire ankle joint versus the competition which relies on a single point of fixation.

- Active Heel Technology
- Accommodates foot from sizes 5 (female) to 14 (male)
- Universal: One boot fits left or right foot
- Easily retrofits your existing Mizuho OSI, Steris®, Maquet® and Schaerer® traction tables
- Eliminate the need for excessive taping
- Patient Weight Capacity: 500 lbs (227 kg)

Choose Your Traction System

Steris OrthoVision® System



#A-93005
Traction Boot



#A-93006
Clamp for Steris

Mizuho OSI Ovation® & Modular Table Systems



#A-93005
Traction Boot



#A-93007
Clamp for Mizuho OSI

Maquet & Schaerer Systems



#A-93000
Traction Boot



#A-93002
Clamp for Maquet
& Schaerer

What's Required



#A-93001
Traction Boot
Disposable
Case of 6



Deluxe Arthroscopic Leg Holder #A-10015

Smooth, one directional locking mechanism enables the caregiver to set rigid fixation with one hand operation. The hinged upper brace accommodates patients while the structurally enhanced design enables the staff to apply pressure both medially and laterally. The ratchet-style securement strap easily adjusts for precise amount of compression. Holds legs up to 35" (89 cm) in circumference. Patient Weight Capacity: 500 lbs (227 kg).

What's Included



#A-10015
Deluxe Arthroscopic
Leg Holder
1 Included



#R-S4133-V
Leg Holder Disposable Pad
3 Included

What's Required



#A-40022 (US)*
Easy Lock® Socket
1 Required

What's Optional



#O-AWLH
Arthroscopic
Well-Leg Holder



#H-009-00-120
Foam Well-Leg Holder



#A-30010 (US)*
Accessory Cart

*Specify Rail Type: US, EU, UK, JP, DEN. See product page for international part numbers.



Classic Arthroscopic Leg Holder #R-734-33

The Arthroscopic Leg Holder immobilizes the thigh and allows for lower leg manipulation during knee arthroscopy. A removable slotted post facilitates the placement of the patient's leg into the holder while the stabilizing bar restricts upper leg movement during knee manipulation. Patient Weight Capacity: 500 lbs (227 kg).

What's Included



#R-734-33
Classic Arthroscopic
Leg Holder
1 Included



#R-S4133-V
Leg Holder Disposable Pad
3 Included

What's Required



#A-40022 (US)*
Easy Lock Socket
1 Required

What's Optional



#O-AWLH
Arthroscopic
Well-Leg Holder



#H-009-00-120
Foam Well-Leg Holder



#R-734-34
Replacement Slide Bar

*Specify Rail Type: US, EU, UK, JP, DEN. See product page for international part numbers.



Leg Holder Disposable Pads #R-S4133-V €€

The Leg Holder Disposable Pad is compatible with all Allen® leg holder systems. These pads utilize two special foam materials which are laminated together. The outer pad is a dense, slightly compressible material which allows for high clamping pressure to stabilize the leg. The inner pad is a soft, fully compressible material for gentle, non-abrasive contact with the skin. The leg pad secures with simple Velcro® fasteners. Length is 37¾" (95.9 cm) and width is 3¾" (8.6 cm). (Sold In Cases Of 12)



Arthroscopic Well-Leg Holder #O-AWLH €€

Easily and safely position the non-affected leg during arthroscopic knee procedures. Help reduce the discomfort and circulatory impairment that can result from the well-leg hanging over the edge of the table. Smooth adjustment and locks securely for precise positioning.

What's Included



#O-AWLH
Arthroscopic
Well-Leg Holder
1 Included



#O-S1007A0.6S
Boot Pad
1 Included

What's Required



#A-40022 (US)*
Easy Lock® Socket
1 Required

What's Optional




#A-30010 (US)*
Accessory Cart

*Specify Rail Type: US, EU, UK, JP, DEN. See product page for international part numbers.



Foam Well-Leg Holder #H-009-00-120

The Foam Well-Leg Holder supports with comfort the non-affected, well-leg during arthroscopic knee procedures, fracture repairs and cast procedures.

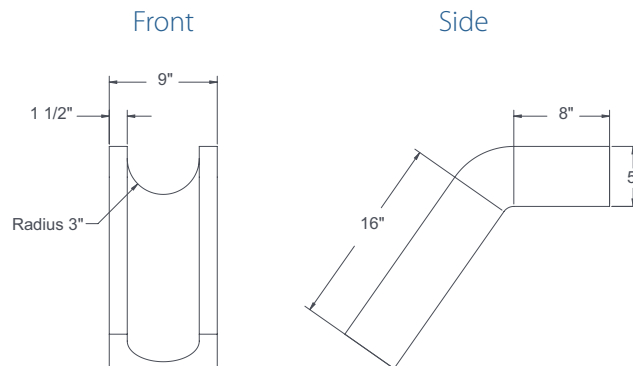

Arthroscopic Stress Post #O-ASP

A levering fulcrum for improved access to the medial compartment of the knee. The 3" x 8" (7.6 cm x 20 cm) post surface is covered with a 1" (2.5 cm) thick, soft pressure management pad for maximum patient safety and comfort.

What's Required



#A-40022 (US)*
Easy Lock® Socket
1 Required

*Specify Rail Type: US, EU, UK, JP, DEN. See product page for international part numbers.



Total Knee Stabilizer #O-TKS

A convenient, stable device to support the flexed leg during total knee procedures. The Total Knee Stabilizer eliminates awkward sandbags and simply mounts to the table rail. The 10" long by 2" diameter (25.5 cm x 5 cm) foot support is covered with a soft, replaceable closed-cell foam pad. The 13" (33 cm) long mounting bar permits vertical height adjustment while the flat 1/4" x 1" (0.6 cm x 2.5 cm) cross section prevents unwanted rotation. Patient Weight Capacity: 500 lbs (227 kg).

What's Required



A-40040 (US)*
Easy Lock® Blade Clamp
1 Required



#O-TKSP
Replacement
Disposable Pad
1 Required

*Specify Rail Type: US, EU, UK, JP, DEN. See product page for international part numbers.



The Prepper Device #O-LPS

Padded cradle comfortably supports patient's leg during prep procedure. Shaft can be rotated to move offset in various positions for convenient access to the leg. Post is 30 1/2" (77.5 cm) long.

What's Required



A-40022 (US)*
Easy Lock Socket
1 Required



#O-LPSP
Replacement
Disposable Pad
1 Required

*Specify Rail Type: US, EU, UK, JP, DEN. See product page for international part numbers.



Radiolucent & Durable



Rectangular



Hourglass



Minor Procedure

Arm & Hand Tables Comparison

Allen® Medical Systems offers an arm and hand table for every need. All of the tables are radiolucent with outstanding durability, backed by a two-year warranty.

- The carbon fiber, Carbon Lights™ tables utilize a frameless design to avoid perimeter interference during radiographic imaging
- The Basic tables utilize a phenolic board with a perimeter frame to allow for radiographic imaging

		#A-30450-C Carbon Lights	#O-AHTR Basic
Rectangular	Weight	5.8 lbs (2.6 kg)	8.8 lbs (4 kg)
	Dimensions	15" x 30 ³ / ₈ " (38 cm x 77 cm)	14" x 32" (35.5 cm x 81 cm)
	Material	Carbon Fiber	Phenolic
	Height Adjustable	Yes	Yes
	Legs	Optional	Integral
	Post Type	Round	Blade
	Clamp Type	Universal Clamp	Rail Clamp
	Pad	Sold Separately	2" (5 cm) Standard Pad

		#A-30650-C Carbon Lights	#O-AHTH Basic
Hourglass	Weight	9.2 lbs (4.2 kg)	11 lbs (5 kg)
	Dimensions	16" x 34" (41 cm x 86 cm)	14" x 31" (35.5 cm x 79 cm)
	Material	Carbon Fiber	Phenolic
	Height Adjustable	Yes	Yes
	Legs	Integral	Integral
	Post Type	Round	Blade
	Clamp Type	Integral (US and EU)	Rail Clamp
	Pad	Sold Separately	2" (5 cm) Standard Pad

		#A-30410-C Carbon Lights	#O-AHMPT Basic
Minor Procedure	Weight	3.4 lbs (1.5 kg)	7.4 lbs (3.4 kg)
	Dimensions	15" x 24" (38 cm x 61 cm)	14" x 26" (35.5 cm x 66 cm)
	Material	Carbon Fiber	Phenolic
	Height Adjustable	Yes	Yes
	Legs	Legless	Legless
	Post Type	Round	Blade
	Clamp Type	Universal Clamp	Rail Clamp
	Pad	Sold Separately	2" (5 cm) Standard Pad



Outer Rail Requires US Rail Clamp for Accessories

Basic Rectangular Table #O-AHTR €€

The Basic Rectangular Table features a telescoping leg that can be easily moved in various positions to provide additional support. The table permits superior access to the surgical site and it attaches to a standard OR table. Includes 1.5" (3.8 cm) pad, for use with 2" (5 cm) OR table pads.

- Weight: 8.8 lbs (4 kg)
- Table Dimensions: 14" x 32" (35.5 cm x 81 cm)
- Patient Weight Capacity: 400 lbs (181 kg)

What's Included



#O-AHTR
Basic Rectangular Table
1 Included



#O-AHTP
1.5" (3.8 cm) Pad
1 Included

What's Required



A-40040 (US)*
Easy Lock® Blade Clamp
2 Required

*Specify Rail Type: US, EU, UK, JP, DEN. See product page for international part numbers.



Table With Leg

Carbon Lights™ Rectangular Table €€

#A-30450-C

The perfect answer for arm and hand procedures. Use without the removable leg when access is paramount. Simply add the optional leg for those procedures requiring additional stability. Attaches via two round 5/8" (1.6 cm) posts. Call Customer Service to inquire about clamps to fit alternate rail sizes. Note: Pad sold separately.

- Carbon Fiber
- Weight: Approx. 5.8 lbs (2.6 kg)
- Table Dimensions: 15" x 30 3/8" (38 cm x 77 cm)
- Patient Weight Capacity: 400 lbs (181 kg)

What's Required *Choose One Pad*



#A-40015 (US Only)
Universal Accessory Clamp
2 Required



#A-30300-Pad-1
1" (2.5 cm) Pad
1 Required

OR



#A-30300-Pad
2" (5 cm) Pad
1 Required

What's Optional



#A-20750
Schlein Hand Positioner
Includes Rail Attachment



#A-30450-LEG
Optional Leg



Outer Rail Requires US Rail Clamp for Accessories

Basic Hourglass Table #O-AHTH

The Basic Hourglass Table features a telescoping leg that can be easily placed in various positions to provide additional support. The table provides additional surface for surgeon's elbows or instrumentation. Attaches to a standard OR table and includes a 1.5" (3.8 cm) pad, for use with 2" (5 cm) OR table pads.

- Weight: 11 lbs (5 kg)
- Table Dimensions: 31 $\frac{3}{8}$ " x 23" *Max* / 14 $\frac{3}{4}$ " *Min* (79.7 x 58.4 *Max* / 37.5 *Min* cm)
- Patient Weight Capacity: 400 lbs (181 kg)

What's Included



#O-AHTH
Basic Hourglass Table
1 Included



#O-AHTHP1.5
1.5" (3.8 cm) Pad
1 Included

What's Required *Choose One*



A-40040 (US)*
Easy Lock® Blade Clamp
2 Required



#O-RC2 (US)*
Rail Clamp
2 Required

*Specify Rail Type: US, EU, UK, JP, DEN, SWISS. See product page for international part numbers.



Carbon Lights™ Hourglass Table #A-30650-C

This full-sized table provides stability and is light enough for easy transport. The adjustable support leg provides additional stability and can be tucked under for easy access. The integral mounting clamps allow vertical adjustment to align with pads on X-ray tops. Note: Pad sold separately.

- Carbon Fiber
- Fits EU & US Rails Only
- Weight: Approx. 9.2 lbs (4.2 kg)
- Table Dimensions: 34" x 16" (26" at end) 86 cm x 41 cm (66 cm at end)
- Patient Weight Capacity: 400 lbs (181 kg)

What's Included



#A-30650-C
Carbon Lights
Hourglass Table
1 Included




#A-30650ER
Elbow Rest
1 Included

What's Required *Choose One*



#A-30502-1-B*
Pad Set for A-30650-C
1" (2.5 cm)
1 Required

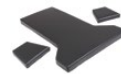

#A-30502-B*
Pad Set for A-30650-C
2" (5 cm)
1 Required



Integral Clamps

What's Optional



#A-30650-Clamp
Replacement Clamp



#A-20750
Schlein Hand Positioner
Includes Rail Attachment



#A-30530-1
Elbow Rest Pads

*Pad sets include table pad and two elbow rest pads



Table With Leg


Carbon Lights™ Minor Procedures Table

#A-30410-C

This small and ultra lightweight table works well for procedures such as carpal tunnel releases and AV fistulas. Has two round 5/8" (1.6 cm) posts. Note: Pad sold separately.


- Carbon Fiber
- Weight: Approx. 3.4 lbs (1.5 kg)
- Table Dimensions: 24" x 15" (61 cm x 38 cm)
- Patient Weight Capacity: 400 lbs (181 kg)

What's Required *Choose One Pad & One Clamp*

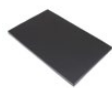

#A-30410-Pad1
1" (2.5 cm) Pad
1 Required

OR

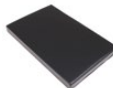

#A-30410-Pad
2" (5 cm) Pad
1 Required



#A-40015 (US Only)
Universal Accessory Clamp
2 Required

OR



#O-RCP2 (US/EU)
Universal Clamp
2 Required

What's Optional



#A-20750
Schlein Hand Positioner
Includes Rail Attachment



#A-30450-LEG
Optional Leg



Outer Rail Requires US Rail
Clamp for Accessories

Basic Carpa® Table #O-AHMPT

#O-AHMPT

The Carpa Table is a legless design with a rigid self-bracing frame. It has a 14" x 26" (35.5 x 66 cm) radiolucent top for maximum access and weighs 7.4 lbs (3.4 kg). The perfect choice for carpal tunnel procedures. Easy drop-in rail attachment. Includes 1.5" (3.8 cm) pad, for use with 2" (5 cm) OR table pads.

- Weight: 7.4 lbs (3.4 kg)
- Table Dimensions: 14" x 26" (35.5 cm x 66 cm)
- Patient Weight Capacity: 400 lbs (181 kg)

What's Included



#O-AHMPT
Basic Carpa Table
1 Included



#O-AHMPT
1.5" (3.8 cm) Pad
1 Included

What's Required



#A-40040 (US)*
Easy Lock® Blade Clamp
2 Required

*Specify Rail Type: US, EU, UK, JP, DEN. See product page for international part numbers.



The Equalizer #R-734-05

Used when three or more digits are placed in traction, The Equalizer distributes the weight of the arm equally between each digit. This stainless steel chain and pulley system is fully autoclavable and comes with an attaching clip. The Equalizer can also be used with the Wrist Fracture System and Weightless Shoulder Suspension System.

Required



Sterile Mesh Finger Trap Disposables
See page 98

What's Optional



#P2000
Thermedx® FluidSmart™
Fluid Management System



#O-LPA
Allen® Easy
Irrigation Tower



#R-734-10
Wrist Arthroscopy
System



Wrist Arthroscopy System #R-734-10

The table mounted Wrist Arthroscopy System is a stable platform and ideal for distraction of the hand, wrist, forearm and elbow. Support extensions mounted off of the tower provide wrist flexion and extension, radial and ulnar deviation and allow for dorsal and volar approaches. Autoclavable.

- Tension adjustments up to 40 lbs (18 kg)
- Allen self-equalizing Finger Traps distribute tension evenly
- Wrist stabilizer provides wrist flexion and extension
- Counter traction device eliminates need for hanging weights

What's Included



#R-734-10
Wrist Arthroscopy System
1 Included

What's Required



#A-40018 (US)*
Simple Clamp
1 Required



Sterile Mesh Finger Trap Disposables
See page 98


What's Optional



#P2000
Thermedx® FluidSmart™
Fluid Management System



#O-LPA
Allen Easy Irrigation Tower

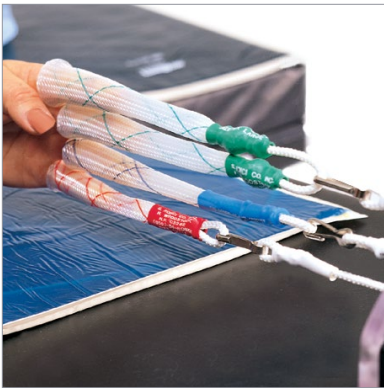

Sterile Mesh Finger Trap Disposables €

Sterile single or double traps. Disposable, for single use only. Five color-coded sizes: Pediatric, Small, Medium, Large, Extra Large.

Sterile Starter Kit

Comes with eight finger traps in four sizes (Extra Large, Large, Medium, Small - two each)

#R-30522-999 Double Starter Kit

#R-30522-888 Single Starter Kit



Pediatric Trap

Eaches - Mix & Match *Minimum Order of 8*

Singles

- #A-30522-8 Single, Extra Large
- #A-30522-2 Single, Large
- #A-30522-4 Single, Medium
- #A-30522-5 Single, Small
- #A-30522-9 Single, Pediatric trap

Doubles

- #A-30522-7 Double, Extra Large
- #A-30522-1 Double, Large
- #A-30522-3 Double, Medium
- #A-30522-6 Double, Small

Cases - Buy In Bulk *8 Per Case*

Singles

- #R-30522-8 Single, Extra Large
- #R-30522-2 Single, Large
- #R-30522-4 Single, Medium
- #R-30522-5 Single, Small
- #R-30522-9 Single, Pediatric Trap

Doubles

- #R-30522-7 Double, Extra Large
- #R-30522-1 Double, Large
- #R-30522-3 Double, Medium
- #R-30522-6 Double, Small



Sizing Kit #R-50522-1 €

Facilitates size selection prior to opening sterile finger traps.



Schlein Hand Positioner #A-20750 €

Supports multiple hand positions. Innovative cleat system can be used with or without weights. Autoclavable.

What's Included



#A-20750
Schlein Hand Positioner
1 Included



Rail Attachment
1 Included

What's Required



Sterile Mesh Finger
Trap Disposables
See page 98

What's Optional



#A-30410-C
Carbon Lights™ Minor
Procedures Table



#A-30450-C
Carbon Lights
Rectangular Table



#A-30650-C
Carbon Lights
Hourglass Table



Malleable Hand Positioner #A-20900

This aluminum malleable hand positioner is perfect for micro nerve reconstruction procedures. As a limited reusable product, the Hand Positioner can withstand more than 20 articulations. It can be sterilized and is radiolucent. (Sold In Cases Of 5)



Portable Orderly Limb Positioner #A-20001

Padded cradle safely supports arm during elbow procedures. Pivoting cradle automatically adjusts to the correct angle required for the selected position.

What's Included



#A-20001
Portable Orderly
Limb Positioner
1 Included



#A-10105
Small Armboard Gel Pad
1 Included



#A-20100
Velcro® Strap
1 Included

Required



#A-40022 (US)*
Easy Lock® Socket
1 Required

*Specify Rail Type: US, EU, UK, JP, DEN. See product page for international part numbers.



Mediq Suomi Oy
PL 115, 02201 Espoo
Puhelin 020 112 1510
asiakaspalvelu@mediq.com, www.mediq.fi